

R3Charge Pre-stained Protein Ladder III, 10-250 kDa

Pack Sizes & Catalogue Numbers

- i. **DS1282120S** – 20 µL, 4 Reactions (Sampling Kit)
- ii. **DS1282120** – 250 µL, 50 Reactions
- iii. **DS1282126** – 6 × 250 µL, 300 Reactions

Product Description

The **R3Charge Pre-stained Protein Ladder III** contains a mixture of **11 prestained proteins** ranging from **10 kDa to 250 kDa** (10, 15, 20, 25, 30, 40, 50, 72, 100, 130, 250 kDa).

- The **72 kDa band** appears **orange**,
- The **25 kDa band** appears **green**,
- All other bands appear **blue**.

This marker is designed for:

- Monitoring protein separation during **SDS-PAGE electrophoresis**,
- Checking **Western blot transfer efficiency** (PVDF, nylon, or NC membranes),
- Estimating the **molecular weight** of proteins.

The ladder is pre-mixed in loading buffer, ready to use. **No heating, dilution, or reducing agent is required before loading.**

Buffer Composition

- 62.5 mM Tris-phosphate buffer (pH 7.5, 25 °C)
- 2% (w/v) SDS
- 5 mM EDTA
- 0.02% (v/v) ProClin300
- 33% (v/v) Glycerol

Usage Instructions


1. Thaw at room temperature and mix thoroughly before use. **Do not boil.**
2. Load **5 μ L per lane** along with experimental samples for SDS-PAGE.
 - For standard 5×1.5 mm wells, 5 μ L is sufficient to visualize clear bands during and after electrophoresis or after transfer.
3. Adjust loading volume (5–10 μ L) according to gel well size and experimental needs.

Storage Conditions

- Store at **–20 °C** for up to **1-year**.
- Store at **4 °C** for up to **3-month** (avoid repeated freeze–thaw cycles).

Notes & Precautions

- Do not run gels excessively long, as bands may diffuse.
 - For frequent use, aliquot after thawing and store at 4 °C for short-term use.
 - Always wear a lab coat and disposable gloves when handling.
-

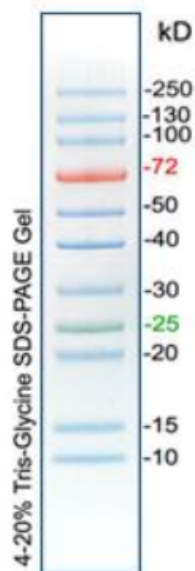

For in vitro research use only, not for diagnostic or therapeutic use. This product is not a medical device.

Contact:

- Support@darwinscience.com